

# Rascal P200



**Mid-range powerchairs**

**Your new route to independence**

# Rascal P200

## Rascal P200 Specification:

Type Class	B / 2
Overall length	109cm / 42.9"
Overall width	73.7cm / 29"
Floor to seat height	57cm (min)
Seat width	45.7cm / 18"
Seat depth	40.6cm / 16"
Maximum carrying capacity	136kgs / 21 stone
Battery voltage and capacity	12V / 55Ah x 2
Maximum speed	6.4km/h / 4mph
Ground clearance	5cm / 2"
Maximum kerb height	5cm / 2"
Range*	up to 26km / 16.2 miles
The maximum safe slope	8°
Turn-around width	1.32m / 4'3"
Colour	Silver

\*Range on full charge and flat ground.

## Main Features:

- Adjustable upholstery on back
- Removable armrests
- Coil spring suspension
- Convenient joystick charging port
- Solid tyres all round
- Flexible seating system
- Swing-away footrests

## Optional Extras:

- Optional accessory elevating legrests
- Optional lighting pack compatible



## Optional Model:

- The P200 Mini is the narrower version with 16" seat width, 26.1" overall frame width and 34Ah batteries


## Rascal P200 Profile:

The Rascal P200 4 mph Rear-Wheel Drive powerchair suitable for indoor and outdoor use comes with 55Ah batteries as standard and features a flexible seating system with a removable seat cushion enabling bespoke adaptations. The capacity to adjust the seat tilt, backrest angle and the horizontal strap tension for total comfort achieves the ideal seating position.

### Coil Suspension



Up to 26km

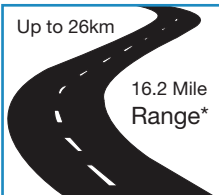


16.2 Mile Range\*

Local Rascal Mobility Retailer:

### Great Turn-around Width

1.30m / 4'3"

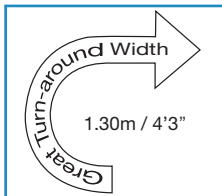

### Swing-away Footrests



Specifications subject to change without prior notice  
E. and O.E- v 2 Final-August 2014