

# Rascal Luxe

Lithium powered scooters




The Lithium Ferro-Phosphate Battery is built to last and is one of the safest Lithium Batteries in the Mobility sector

**RASCAL**  
**EXCLUSIVE**

# Rascal Luxe Profile:

The **Rascal Luxe** 4mph **Lithium Ferro-Phosphate** powered transportable scooter is sleek and durable. With a lightweight aluminium frame, designed for ergonomic “All day” comfort offering increased legroom that accommodates a range of heights from 4’11” - 6’4” plus.

This premium scooter features a supremely comfortable deluxe seat, encouraging **100°** posture perfect alignment of the feet, legs and back position for a supportive journey. The space around the tiller, offers significant leg and knee room that promotes even weight distribution and enhances the stability alongside the adjustable all-round suspension.

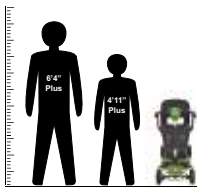

The **Rascal Luxe** benefits from a quality design that showcases styling, an impressive ride up to 21.4 stone user weight and reliable range up to 21.7 miles (90kg user weight). Luxuriously developed comfort and manoeuvrability for every journey with this class leading scooter!

## Main Features:

Two easy access charging points on tiller & battery box

Range 21.7 miles / 34.9km (90kg user weight)\* with **27Ah** battery & Rapid 5A charger

360° full beam light switch

USB port

360° Postural Seating

Fold down backrest with handy pouch

Sleek lightweight, aluminium frame

Advanced stability with impressive legroom

Ergonomic “All day” comfort

Efficient 360° LED lighting

All-round Suspension with exceptional ride comfort

8” Dual compound puncture proof tyres

Mass weight of 27Ah Battery Pack 6.7kg / 14.7lbs

**RASCAL Life** **BRITISH MADE**  
Medical / Military grade **Lithium Ferro-Phosphate Battery**

## Where exceptional design meets safety...



"See and be seen" with the **360° LED lighting**, highly visible from any angle



Useful **tiller lock** helps with lifting, **weighs just 10.8kg** and super portable for easy storage



IPX5\* waterproof rating.  
Adjustable **all-round suspension**

\* Ingress Protection Rating



Practical **USB port**



The **Rascal Luxe** with high quality **Lithium Ferro-Phosphate (LiFePO<sub>4</sub>) Battery**.

Relax and enjoy the ride in the knowledge that the Lithium Ferro-Phosphate (LiFePO<sub>4</sub>) Battery pack has an increased efficiency and performance to stand up to everyday life.

Plus the **LiFetime Battery Guarantee**† for peace of mind!

[Scan to visit website](#)



## Truly enhance the experience and freedom to travel

### Explorers' Choice - *Have the Best of Both Worlds*



Purchase the 11.6Ah battery option & 2A charger<sup>^</sup>

Battery weighs just 4.1kg / 9lbs

Range up to 10.3 miles / 16.5km (90kg user weight)\*

#### Air travel & CAA Lithium Battery Legislation

Airline Compliant 11.6Ah battery meets I.A.T.A. UN38.3 Passenger Aircraft Transportation requirements for overseas travel and can be taken on a plane. (Subject to Airline Regulation and Restrictions)

† LiFetime Battery Guarantee T&Cs apply

<sup>^</sup> Speed charge & save time - charge between 10-90% in less than 2 hours with a 5A charger

# Rascal Luxe

**LiFetime Battery  
Guarantee\***

**Specification:**      **Luxe 27Ah**      **Luxe 11.6Ah**  
with Rapid 5A Charger      with 2A Charger

Type Class	B / 2	✓
Overall length	103.5cm / 40.7"	✓
Overall width	49.5cm / 19.5"	✓
Overall height	94cm / 37"	✓
Floor to seat height	46-50cm (min-max)	✓
Seat width	44cm / 17.3"	✓
Seat depth	41cm / 16.1"	✓
Folded Dimensions LxWxH	103.5cm / 49.5cm / 71cm (40.7" / 19.5" / 27.9")	✓
Maximum carrying capacity	136kgs / 21.4 stone	✓
Battery voltage and capacity	25.6V / 27Ah	25.6V / 5.8Ah x 2
Maximum speed	6.4km/h / 4mph	✓
Ground clearance	4.5cm / 1.8"	✓
Maximum kerb height	5cm / 2"	✓
Range* (Max kg user weight)	up to 25.8km / 16 miles	up to 12.2km / 7.6 miles
Range* (90kg user weight)	up to 34.9km / 21.7 miles	up to 16.5km / 10.3 miles
The maximum safe slope	6°	✓
Turn-around width	1.25m / 4'1"	✓
Colour	Champagne / Jade Green / Royal Blue	
Single Heaviest component (front section)	10.8kgs / 23.8lbs	✓
Product weight (Excluding battery)	32kgs / 70.5 lbs	✓
(Including battery)	38.7kgs / 85.3 lbs	36.1kgs / 79.6 lbs

\*Range on full charge and flat ground

\*See full Tech Page on Website & Owner's Manual for details

## Practical Features:



† LiFetime Battery Guarantee T&Cs apply

**ELECTRIC  
MOBILITY**



Specifications subject to change without prior notice  
Images are for illustrative purposes only  
E. and O. E- v4 April 2026

electricmobility.co.uk